



Care@Home™

C8000 Series Control Panel

Getting Started


For Cellular,
Ethernet,
Wi-Fi and
BLE



The Best Location for the C8000



Good cellular
reception

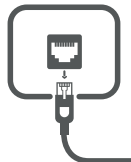

Near network
port



Convenient for
resident to access



Where speakerphone can be
heard throughout premises


Setting Up the C8000



1 Place thumb in oval hole and lift off the cover



3 Insert the backup battery

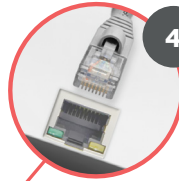

4 Plug in Ethernet cable and attach to network port



2 Insert the SIM




6 Turn the power switch on




5 Close the back cover and secure it with the provided screw. Connect the power cable to the power source.

Using the C8000

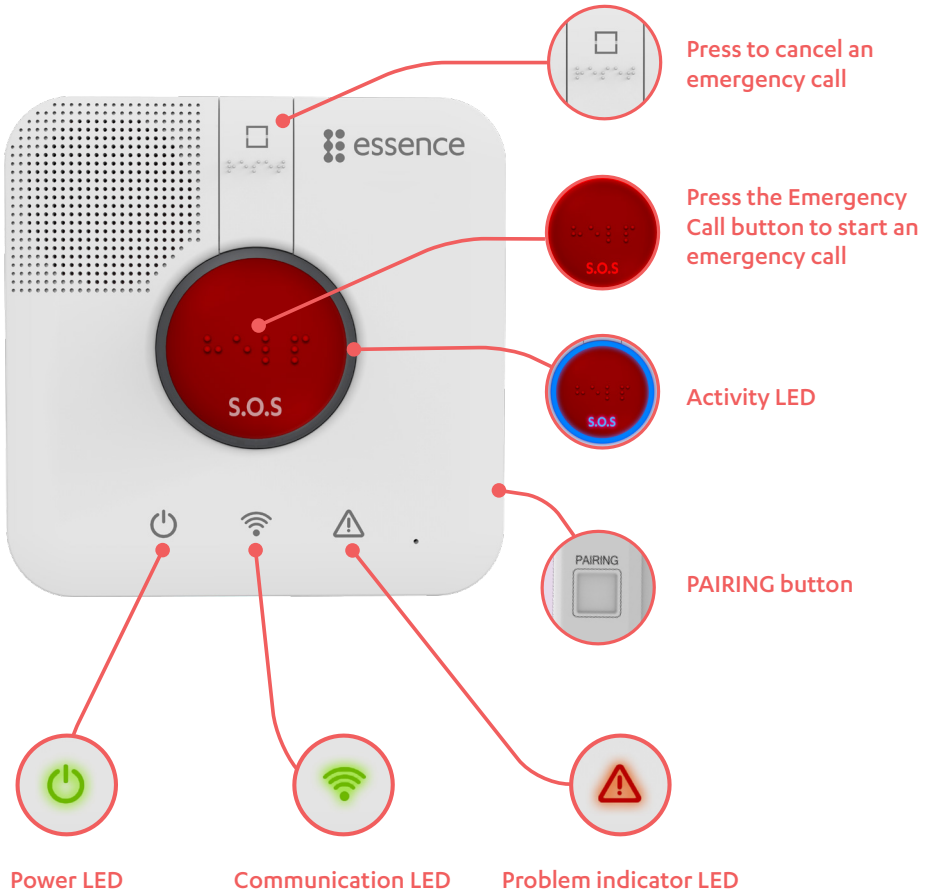

The Emergency Pendant in Action

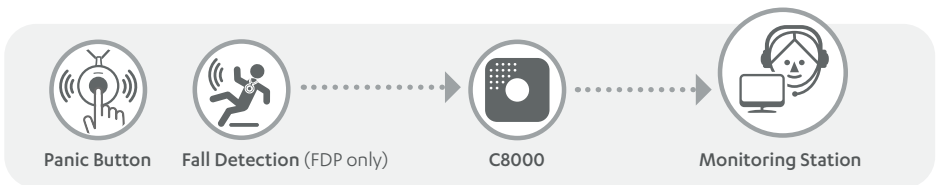

The EPP/FDP blinks **red** when the panic button is pressed
The FDP blinks **red** when a fall is detected

Emergency Pendant Options

The C8000 is also available with the following pendant options:

Emergency Pendant Plus (EPP) without fall detection



Fall Detection Pendant (FDP) including fall detection



Setting Up the Emergency Pendant

1


Insert the battery into the EPP/FDP*>



2

Press the PAIRING button for 3 seconds


The Activity LED ring blinks blue

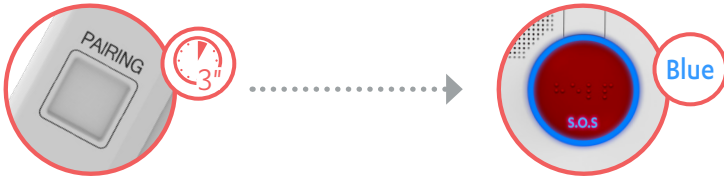

3

Press the PANIC button

"Device added successfully"

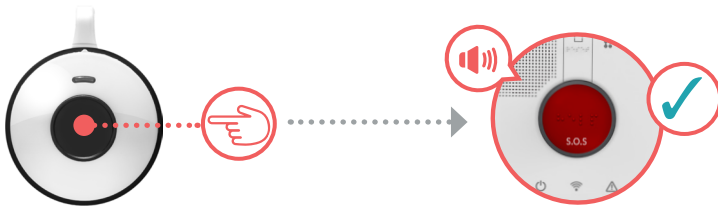

*In some kit configurations, the EPP/FDP is supplied with the battery already installed

